

TECHNICAL SPECIFICATIONS

STENT SPECIFICATIONS			DELIVERY SYSTEM	
Drug/Polymer		Sirolimus/PLGA	Catheter type	Rapid exchange
Drug load		1.4 µg/mm²	Guiding catheter compatibility	5F
Material		Cobalt chromium L605	Guidewire compatibility	0.014"
Strut thickness	Ø 2.25-4.0 mm	0.0026" (65 µm)	Proximal shaft	Stainless steel, PTFE coating
	Ø 4.50-5.0 mm	0.0042" (107 µm)	Coat distal shaft	Hydrophilic coating
Side branch access (Ø 4.0 mm)		5.3 mm	Proximal shaft size	1.9F
Radial strength		28psi	Distal shaft size	2.8F (Ø 2.25 – 4.00 mm)
Foreshortening		< 2.5%		3.0F (Ø 4.50 & 5.00 mm)
Mechanical recoil		< 6%	Nominal pressure	8atm
			RBP	19atm (Ø 2.25mm); 18atm (Ø 2.5-2.75mm); 17atm (Ø 3.0-3.25mm); 16atm (Ø 3.5-4.0mm); 15tam (Ø 4.5-5.0mm)

COMPLIANCE CHART

Pressure (atm)	Balloon Diameter (mm)								
	2.25	2.50	2.75	3.00	3.25	3.50	4.00	4.50	5.00
8 (NOM)*	2.25	2.50	2.75	3.00	3.25	3.50	4.00	4.50	5.00
9	2.27	2.53	2.78	3.03	3.28	3.55	4.04	4.56	5.05
10	2.29	2.56	2.81	3.06	3.32	3.59	4.09	4.60	5.09
11	2.32	2.59	2.84	3.09	3.35	3.62	4.12	4.64	5.13
12	2.34	2.62	2.86	3.12	3.38	3.65	4.15	4.68	5.16
13	2.36	2.65	2.89	3.14	3.41	3.69	4.19	4.72	5.20
14	2.38	2.68	2.91	3.17	3.44	3.72	4.23	4.76	5.23
15 (RBP)**	2.40	2.71	2.94	3.20	3.47	3.76	4.27	4.81	5.27
16 (RBP)**	2.42	2.74	2.97	3.22	3.50	3.80	4.31	-	-
17 (RBP)**	2.45	2.78	3.00	3.25	3.53	-	-	-	-
18 (RBP)**	2.47	2.82	3.03	-	-	-	-	-	-
19 (RBP)**	2.49	-	-	-	-	-	-	-	-

* Nominal Pressure. The nominal in vitro device specifications do not take into account any lesion resistance. ** Rated Burst Pressure. Do not exceed RBP

ORDERING INFORMATION

Stent Ø (mm)	Stent Length (mm)										
	8	10	13	16	18	23	28	33	38	43	48
2.25	ELM22508	ELM22510	ELM22513	ELM22516	ELM22518	ELM22523	ELM22528	ELM22533	ELM22538	ELM22543	ELM22548
2.50	ELM25008	ELM25010	ELM25013	ELM25016	ELM25018	ELM25023	ELM25028	ELM25033	ELM25038	ELM25043	ELM25048
2.75	ELM27508	ELM27510	ELM27513	ELM27516	ELM27518	ELM27523	ELM27528	ELM27533	ELM27538	ELM27543	ELM27548
3.00	ELM30008	ELM30010	ELM30013	ELM30016	ELM30018	ELM30023	ELM30028	ELM30033	ELM30038	ELM30043	ELM30048
3.25	ELM32508	ELM32510	ELM32513	ELM32516	ELM32518	ELM32523	ELM32528	ELM32533	ELM32538	ELM32543	ELM32548
3.50	ELM35008	ELM35010	ELM35013	ELM35016	ELM35018	ELM35023	ELM35028	ELM35033	ELM35038	ELM35043	ELM35048
4.00	ELM40008	ELM40010	ELM40013	ELM40016	ELM40018	ELM40023	ELM40028	ELM40033	ELM40038	ELM40043	ELM40048
4.50	-	-	ELM45013	ELM45016	ELM45018	ELM45023	ELM45028	ELM45033	ELM45038	-	-
5.00	-	-	ELM50013	ELM50016	ELM50018	ELM50023	ELM50028	ELM50033	ELM50038	-	-

For more information please visit our website at <https://orbusneich.com> or contact us:

Corporate Headquarters

Unit 303 & 305, 3/F, Building 20E
Hong Kong Science Park
Shatin, N.T., Hong Kong, China
Phone +852.2802.2288
Fax +852.2507.3532

EMEA Regional Headquarters

Drs. W. van Royenstraat 5
3871 AN Hoevelaken
The Netherlands
Phone +31.33.254.1150
Toll Free Phone 00800.0254.1150*


* Only for Belgium, Denmark, France, Germany, Ireland, Netherlands, Norway, Sweden and UK

eucaLimus™ is manufactured in Germany.

©2024 OrbusNeich Medical Group Holdings Limited or its affiliates. All rights reserved.
Not available for sale in the USA.

G-70-2387 Rev03

eucaLimus and OrbusNeich are trademarks of OrbusNeich Medical Group Holdings Limited or its affiliates.
All data and photos on file. Illustrations are not to be considered as engineering drawings or photos.

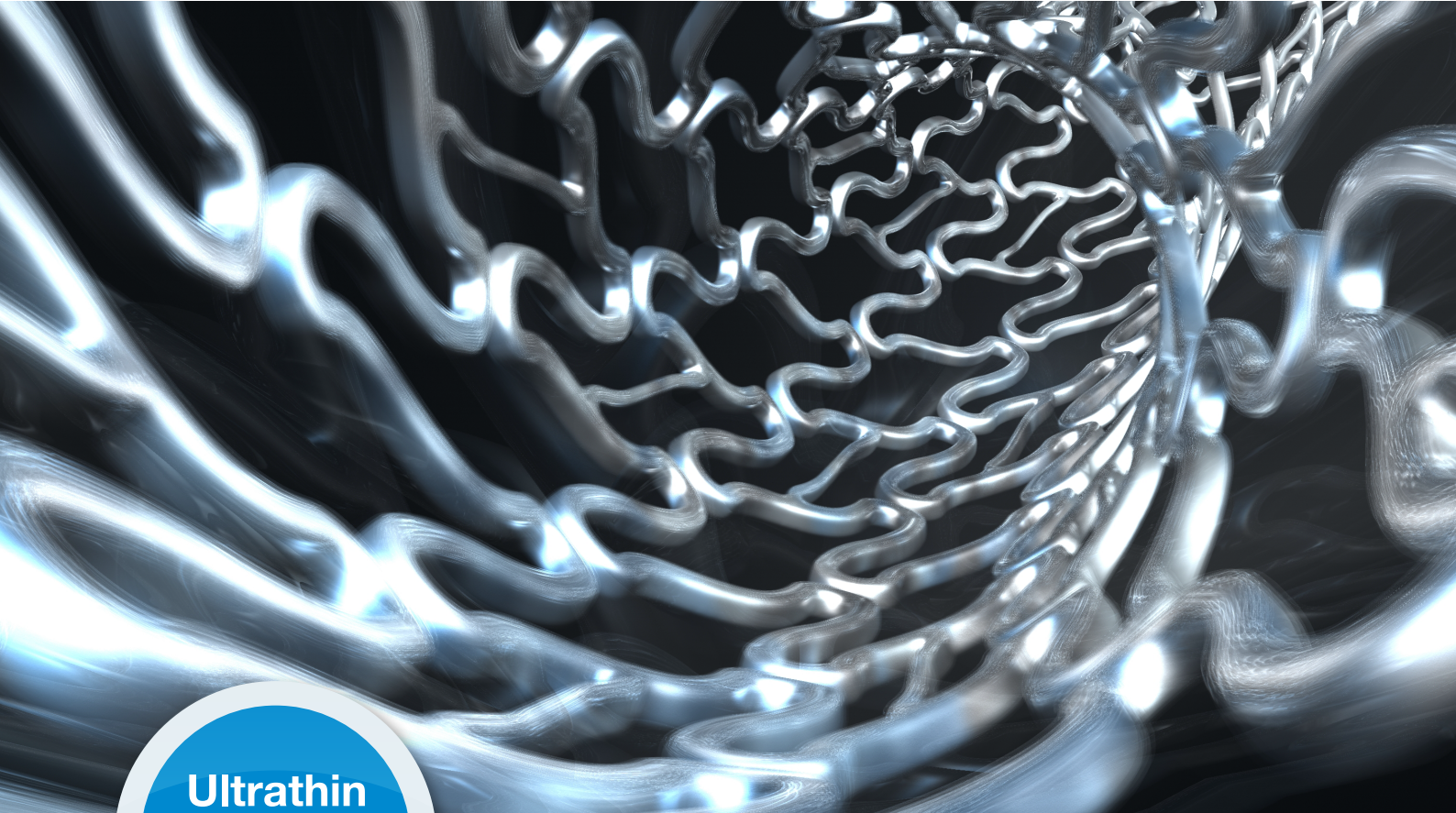


Ultrathin
Strut
65µm


Excellent Deliverability,
Optimal Therapy and Support

Drug/Polymer

Sirolimus

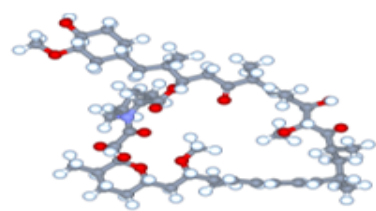


PLGA 85/15 polylactide-co-glycolide

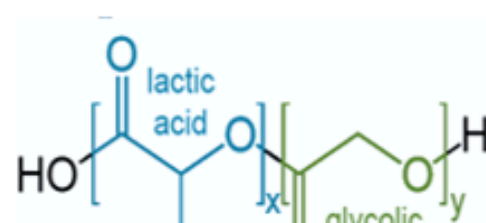

- Sirolimus drug for inhibition of smooth muscle proliferation and intimal hyperplasia
- Drug load is 1.4 µg/mm²
- Drug is eluted via PLGA polymer that degrades within 10-13 weeks period


Coating Technology

Coating Integrity

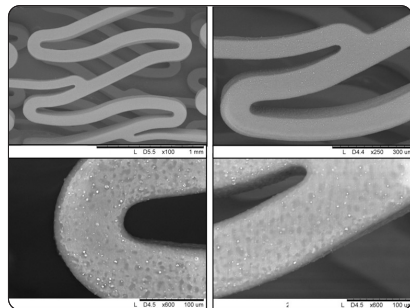


Coating remains intact at maximum post dilatation expansion of 4.5 mm

Asymmetrical Coating

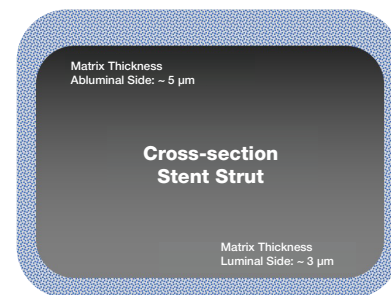


Asymmetrical coating facilitates enhanced efficiency in drug delivery

Ultrathin Strut of 65 µm

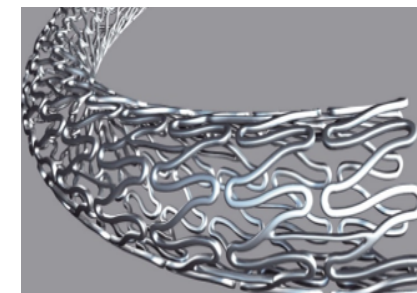


65 µm Strut Thickness

- Excellent deliverability
- Low crossing profile
- Easier crossability

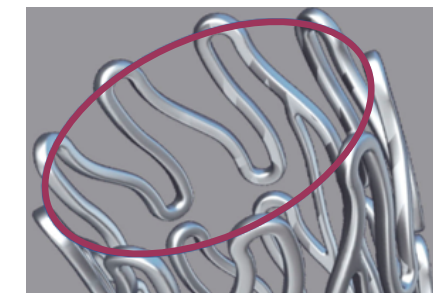


9 Crown Structure

- Good vessel coverage
- Good radial force
- Low recoil & foreshortening

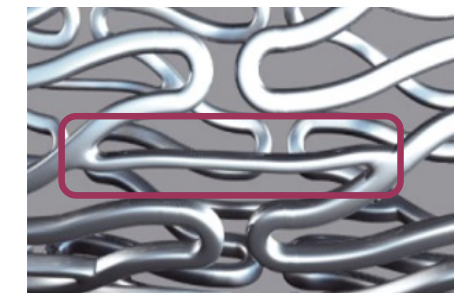


Interlinks

Thin flexible interlinks for enhanced vessel conformability

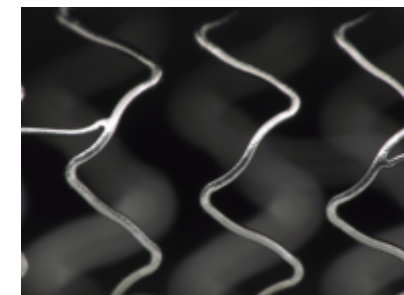

Side Branch Access

Large cell opening facilitates easy device access to side branches

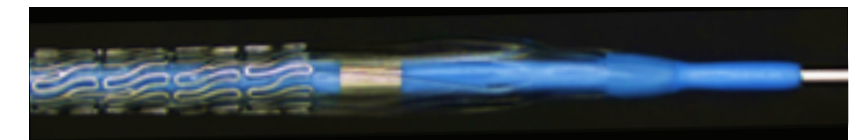

Crossing Profile


Low crossing profile for reachability in complex anatomies

Wide Range of Sizes

- Diameter from 2.25 mm up to 5.0 mm
- Stent length from 8 mm up to 48 mm


Drug Release / Polymer Degradation

Sirolimus Elution Time for Eucalimus

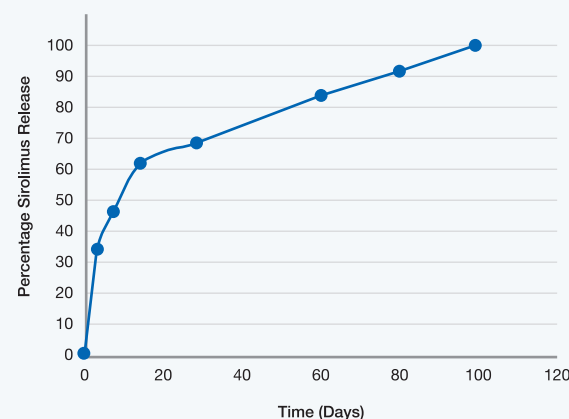


Data on file.

Drug Release and Polymer Degradation

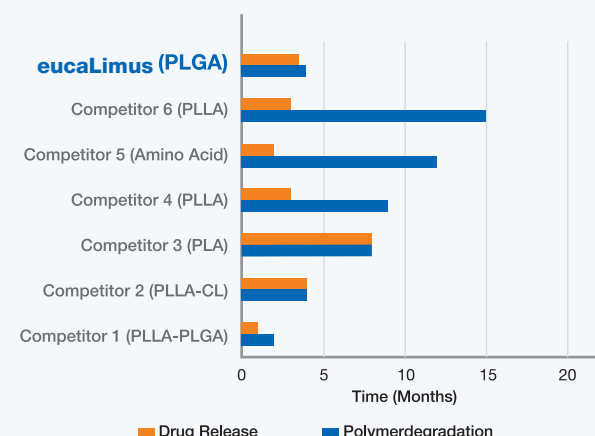

Data on file.

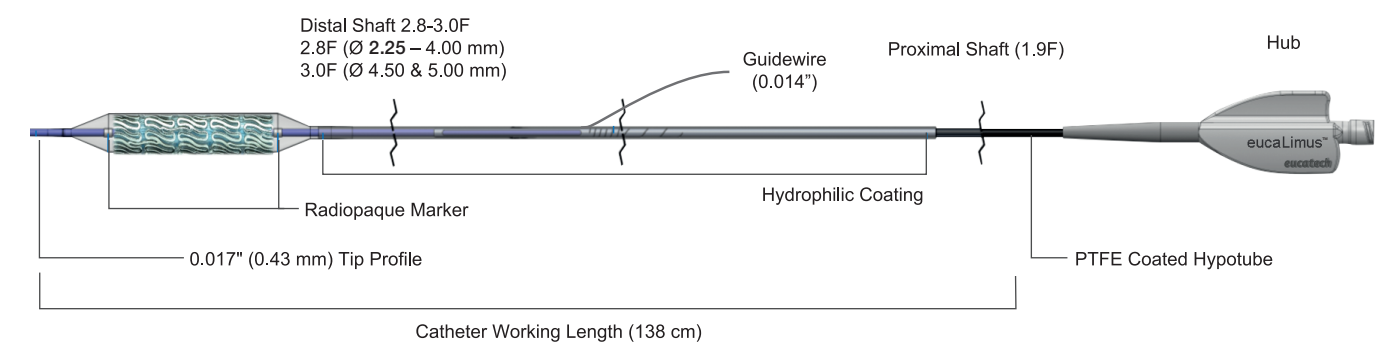

Stent Delivery System

- Enhanced balloon delivery system for pushability and trackability
- High performance hydrophilic coating to minimize friction
- Laser cut hypotube 1:1 transmission zone from hypotube to distal part
- Dual tapered atraumatic tip

36 MONTHS SHELF LIFE

e-beam Sterilization

- Dry sterilization
- Shelf life of 3 years