

Ordering Information


Non-EU Manufacturing Catalogue Numbers

Balloon Diameter (mm)	Balloon Working Length (mm)				
	8	10	12	15	18
1.75	117-1708J	117-1710J	117-1712J	117-1715J	117-1718J
2.0	117-2008J	117-2010J	117-2012J	117-2015J	117-2018J
2.25	117-2208J	117-2210J	117-2212J	117-2215J	117-2218J
2.5	117-2508J	117-2510J	117-2512J	117-2515J	117-2518J
2.75	117-2708J	117-2710J	117-2712J	117-2715J	117-2718J
3.0	117-3008J	117-3010J	117-3012J	117-3015J	117-3018J
3.25	117-3208J	117-3210J	117-3212J	117-3215J	117-3218J
3.5	117-3508J	117-3510J	117-3512J	117-3515J	117-3518J
3.75	117-3708J	117-3710J	117-3712J	117-3715J	117-3718J
4.0	117-4008J	117-4010J	117-4012J	117-4015J	117-4018J
4.5	117-4508J	117-4510J	117-4512J	117-4515J	117-4518J
5.0	117-5008J	117-5010J	117-5012J	117-5015J	117-5018J

Compliance Chart

Pressure atm (kPa)	Balloon Diameter (mm)												
	1.75	2.0	2.25	2.5	2.75	3.0	3.25	3.5	3.75	4.0	4.5	5.0	
6 (608)	1.68	1.90	2.14	2.40	2.61	2.86	3.10	3.33	3.54	3.79	4.26	4.70	
8 (811)	1.70	1.93	2.17	2.44	2.66	2.91	3.15	3.39	3.61	3.86	4.34	4.80	
10 (1013)	1.73	1.97	2.21	2.47	2.70	2.95	3.20	3.44	3.68	3.93	4.42	4.90	
NOM* 12 (1216)	1.75	2.00	2.25	2.50	2.75	3.00	3.25	3.50	3.75	4.00	4.50	5.00	
14 (1419)	1.77	2.03	2.29	2.53	2.80	3.05	3.30	3.56	3.82	4.07	4.58	5.10	
16 (1621)	1.80	2.07	2.33	2.56	2.84	3.09	3.35	3.61	3.89	4.14	4.66	5.20	
RBP** 18 (1824)	1.82	2.10	2.36	2.60	2.89	3.14	3.40	3.67	3.96	4.21	4.74	5.30	
RBP** 20 (2027)	1.84	2.13	2.40	2.63	2.93	3.18	3.45	3.72	4.03	4.27	4.82	5.40	
22 (2229)	1.87	2.16	2.44	2.66	2.98	3.23	3.50	3.78	4.10	4.34	4.90	5.49	
24 (2432)	1.89	2.20	2.48	2.69	3.02	3.27	3.55	3.83	4.17	4.41	4.98	5.59	
26 (2634)	1.92	2.23	2.51	2.72	3.07	3.32	3.59	3.89	4.24	4.48	5.06	5.69	

* Nominal Pressure. The nominal in vitro device specifications do not take into account any lesion resistance.
** Rated Burst Pressure. Do not exceed RBP

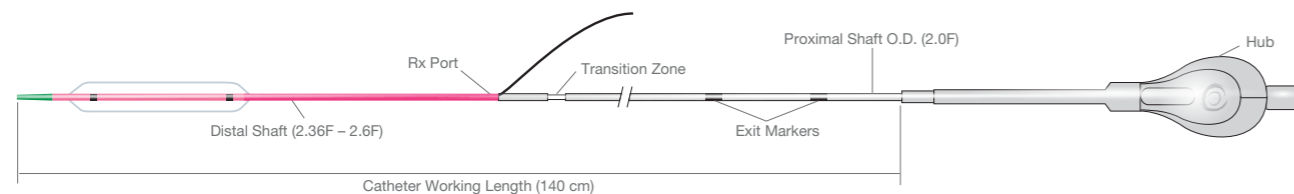

For more information please visit our website at www.OrbusNeich.com or contact us:

Corporate Headquarters

Unit 303 & 305, 3/F, Building 20E
Hong Kong Science Park
Shatin, N.T., Hong Kong, China
Phone +852.2802.2288
Fax +852.2507.3532

EMEA Regional Headquarters

Drs. W. van Royenstraat 5
3871 AN Hoevelaken
The Netherlands
Phone +31.33.254.1150 Toll Free Phone 00800.0254.1150*

* Only for Belgium, Denmark, France, Germany, Ireland, Netherlands, Norway, Sweden and UK

©2020 OrbusNeich Medical Company Limited or its affiliates. All rights reserved.

Not available for sale in the USA.

G-70-0201 Rev04

Sapphire[®] and OrbusNeich are registered trademarks of OrbusNeich Medical Company Limited or its affiliates. All data and photos on file. Illustrations are not to be considered as engineering drawings or photos.




True Controlled Compliance and
Ultimate Crossability for a
Wide Spectrum of Vessels


Ultimate Crossability and Tip Entry Profile

Sapphire II NC Coronary Dilatation Catheter is a well-balanced and truly non-compliant balloon catheter engineered to cross tight and calcified lesions.

Crossing Profile Comparison

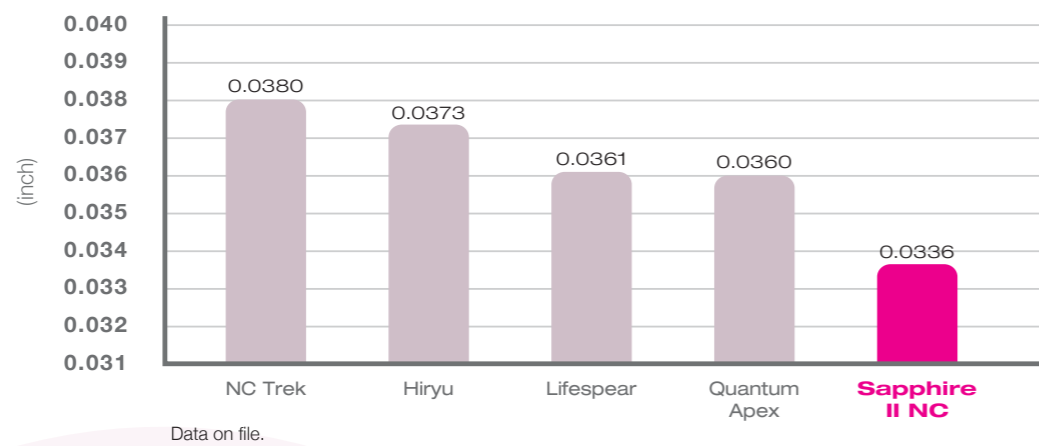


Enhanced distal Z-Tip with advanced laser welding technology for “ZERO” transition with conventional steerable guidewires

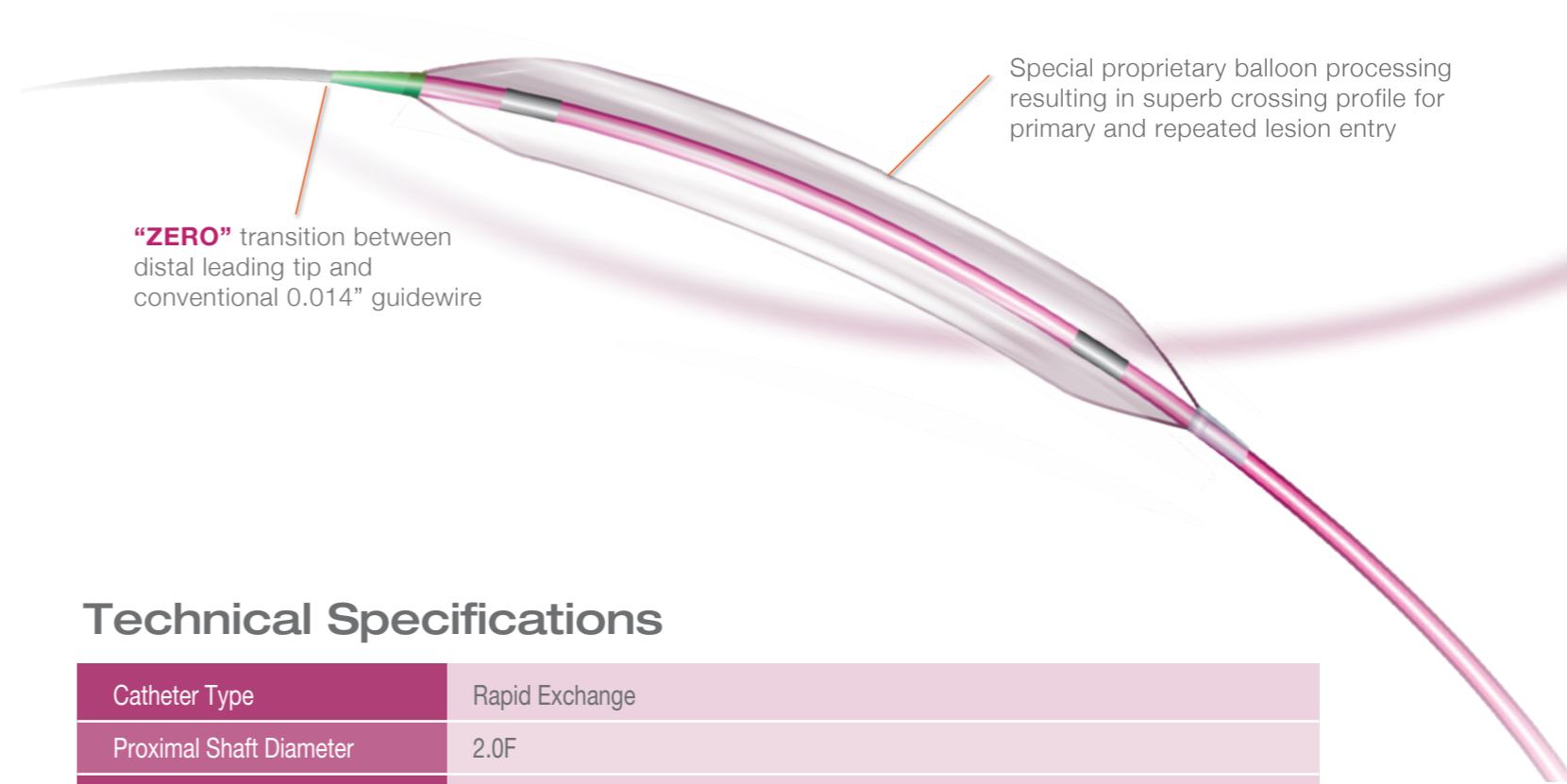

True Controlled Compliance

True non-compliance and minimal **overstretch for reliability** and size accuracy. Sapphire NC **brings** advanced material **engineering for** optimum reliability and deliverability.

Compliance Comparison from Nominal Pressure to RBP

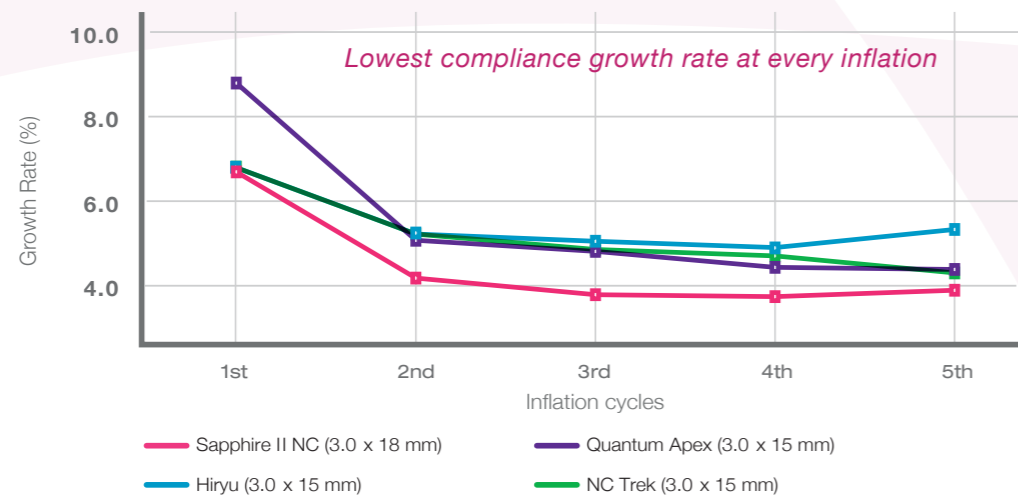

Technical Specifications

Catheter Type	Rapid Exchange
Proximal Shaft Diameter	2.0F
Distal Shaft Diameter	2.36F (Ø 1.75 mm – 2.0 mm); 2.55F (Ø 2.25 mm – 3.5 mm); 2.6F (Ø 3.75 mm – 5.0 mm)
Catheter Working Length	140 cm
Leading Tip Length	1.5 mm (Ø 1.75 mm); 2.0 mm (Ø 2.0 mm – 3.0 mm); 2.5 mm (Ø 3.25 mm – 5.0mm)
Balloon Material	Nylon
Balloon Folds	3 (Ø 1.75 mm – 3.0 mm); 5 (Ø 3.25 mm – 5.0 mm)
Compliance	Controlled-compliant
Crossing Profile	0.0336" (Ø 3.0 mm)
Marker Bands	2
Coating	Hydrophilic (distal tip to guidewire exit port), Hydrophobic (guidewire lumen)
Nominal Pressure	12 atm
Rated Burst Pressure	20 atm (Ø 1.75 mm – 4.0 mm); 18 atm (Ø 4.5 mm – 5.0 mm)
Guiding Catheter Compatibility	5F (1.75 mm – 4.0 mm), 6F (4.5 mm – 5.0 mm)