

Ordering Information (EU Manufacturing Catalogue Numbers)

Ø (mm)	Length (mm)		
	10	15	20
1.75	631751010	631751510	631752010
2.00	632001010	632001510	632002010
2.25	632251010	632251510	632252010
2.50	632501010	632501510	632502010
2.75	632751010	632751510	632752010
3.00	633001010	633001510	633002010
3.50	633501010	633501510	633502010
4.00	634001010	634001510	634002010


Compliance Chart

Pressure (atm)	Balloon Diameter (mm)							
	1.75	2.00	2.25	2.50	2.75	3.00	3.50	4.00
2	1.64	1.87	2.13	2.35	2.55	2.79	3.20	3.67
4	1.66	1.90	2.15	2.38	2.59	2.83	3.26	3.74
6	1.68	1.92	2.18	2.41	2.63	2.87	3.32	3.80
8	1.71	1.95	2.20	2.44	2.67	2.92	3.38	3.87
10	1.73	1.97	2.23	2.47	2.71	2.96	3.44	3.93
12 (NOM)*	1.75	2.00	2.25	2.50	2.75	3.00	3.50	4.00
14	1.77	2.03	2.27	2.53	2.79	3.04	3.56	4.07
16	1.79	2.05	2.30	2.56	2.83	3.08	3.62	4.13
18	1.82	2.08	2.32	2.59	2.87	3.13	3.68	4.20
20 (RBP)**	1.84	2.10	2.35	2.62	2.91	3.17	3.74	4.26
22	1.86	2.13	2.37	2.65	2.95	3.21	3.80	4.33
24	1.88	2.15	2.40	2.68	2.99	3.25	3.86	4.39
26	1.90	2.18	2.42	2.71	3.03	3.29	3.92	4.46

* Nominal Pressure. The nominal in vitro device specifications do not take into account any lesion resistance. ** Rated Burst Pressure. Do not exceed RBP.

Technical Specifications

Scoreflex® TRIO	
Catheter type	Rapid exchange
Guidewire lumen diameter	0.014"
Balloon material	Nylon blend
Balloon compliance	Non-compliant
Coating (distal shaft and tip)	Hydro-X hydrophilic coating
Coating (guidewire lumen)	Invio hydrophobic coating
Crossing profile (Ø3.0 mm)	0.0334" / 0.86 mm
Nominal pressure (NOM)	12 atm
Rated burst pressure (RBP)	20 atm

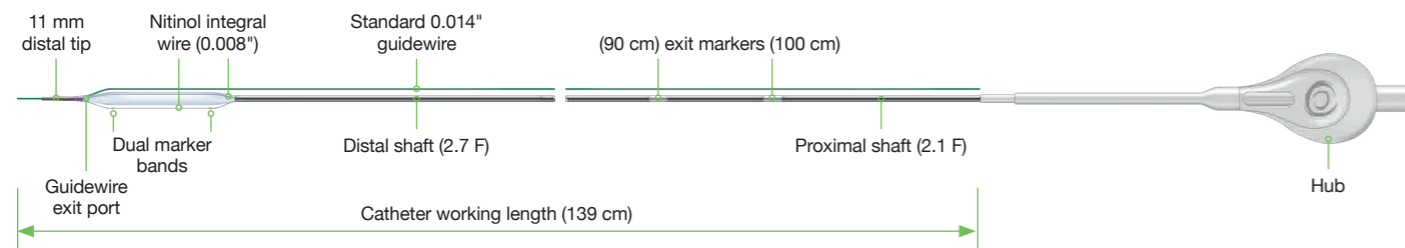

For more information please visit our website at <https://orbusneich.com/> or contact us:

Corporate Headquarters
Units 303 & 305, 3/F, Building 20E
Hong Kong Science Park
Shatin, N.T., Hong Kong, China
Phone +852.2802.2288
Fax +852.2507.3532

EMEA Regional Headquarters
Drs. W. van Royenstraat 5
3871 AN Hoevelaken
The Netherlands
Phone +31.33.254.1150
Toll Free Phone 00800.0254.1150*

*Only for Belgium, Denmark, France, Germany, Ireland, Netherlands, Norway, Sweden and UK

©2025 OrbusNeich Medical Group Holdings Limited or its affiliates. All rights reserved. Scoreflex® Trio is not approved for sale in the USA.

G-70-2476 Rev01

OrbusNeich and Scoreflex are registered trademarks of OrbusNeich Medical Group Holdings Limited or its affiliates. All data and photos on file. Illustrations are not to be considered as engineering drawings or photos.




Greater Scoring Depth, Safety and Deliverability

scoreflex® TRIO

Tri-Scoring Innovation: Greater Scoring Depth, Safety and Deliverability Aligned



Tri-scoring Greater Scoring Depth¹

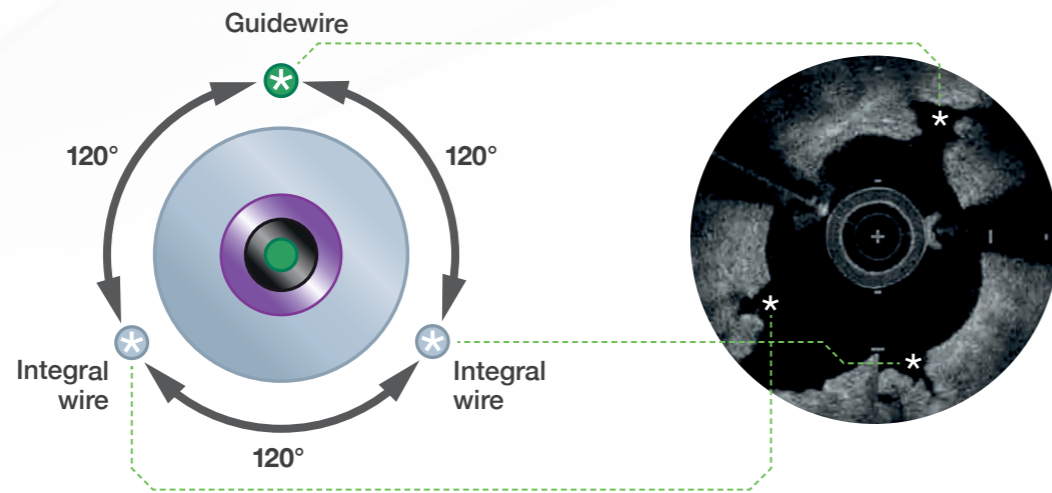


3 distinctive points of plaque modification




Non-compliant Minimal Dog-bone Effect

Dog-bone effect at RBP*

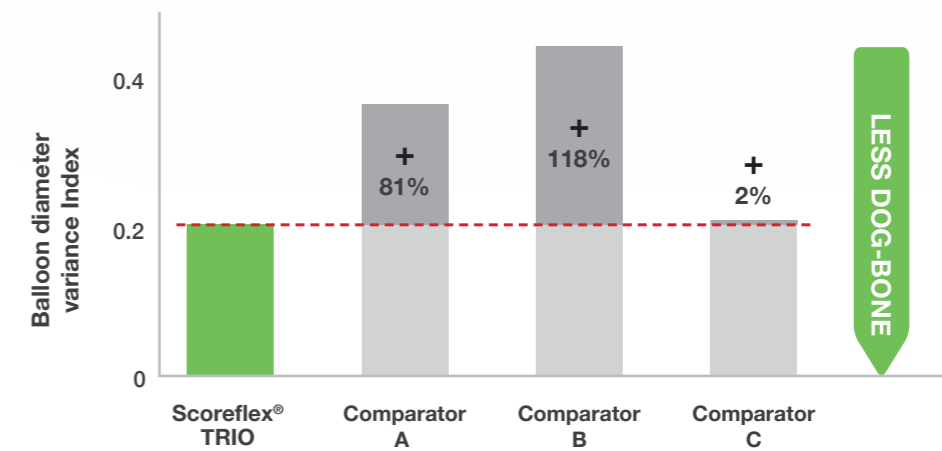


Deliverability

Competitive Crossing Profile*

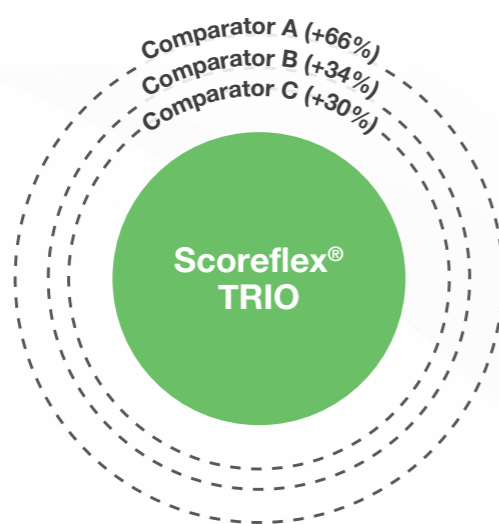


Competitive crossing profile*

Balloon Flexibility*

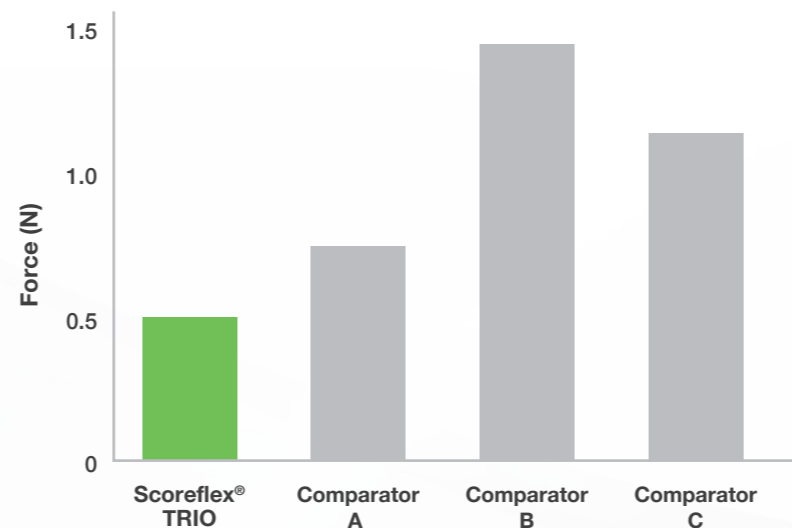

TRIO requires less force for negotiating tortuous anatomies

Reachability²

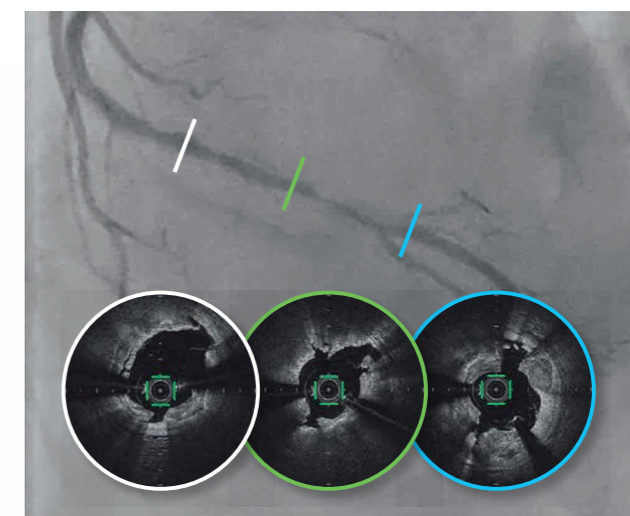

Effective plaque modification seen at distal vessel

* Data on file

¹ Background OCT image courtesy of Dr. Yasuaki Koga, Oita Prefectural Hospital, Japan

² RCA angiographic and OCT images courtesy of Drs. Takanobu Mitarai and Yuki Ishibashi, St. Marianna University Hospital, Japan