

Ordering Information

SZ Manufacturing Catalogue Numbers

Balloon Diameter (mm)	Balloon Working Length (mm)						New Sizes
	15	20	40	80	120	180	240
1.5	L540-15015	L540-15020	L540-15040	L540-15080	L540-15120	L540-15180	L540-15240
2.0	L540-20015	L540-20020	L540-20040	L540-20080	L540-20120	L540-20180	L540-20240
2.5	-	L540-25020	L540-25040	L540-25080	L540-25120	L540-25180	L540-25240
3.0	L540-30015	L540-30020	L540-30040	L540-30080	L540-30120	L540-30180	L540-30240
3.5	-	L540-35020	L540-35040	L540-35080	L540-35120	L540-35180	L540-35240
4.0	L540-40015	L540-40020	L540-40040	L540-40080	L540-40120R*	L540-40180	L540-40240
4.5	-	L540-45020	L540-45040	L540-45080	L540-45120R*	-	-
5.0	L540-50015	L540-50020	L540-50040	L540-50080	L540-50120R*	L540-50180	L540-50240
5.5	-	L540-55020	L540-55040	L540-55080	L540-55120R*	-	-
6.0	L540-60015	L540-60020	L540-60040	L540-60080	L540-60120R*	L540-60180	L540-60240


*R denotes sizes using a modified balloon material subject to domestic approval of each territory

Compliance Chart

Pressure atm	Balloon Diameter (mm)									
	1.5	2.0	2.5	3.0	3.5	4.0	4.5	5.0	5.5	6.0
2	1.41	1.88	2.37	2.82	3.27	3.70	4.22	4.67	5.09	5.60
4	1.43	1.90	2.40	2.86	3.32	3.76	4.27	4.74	5.17	5.68
6	1.45	1.93	2.42	2.89	3.36	3.82	4.33	4.80	5.26	5.76
8	1.47	1.95	2.45	2.93	3.41	3.88	4.39	4.87	5.34	5.84
10	1.48	1.98	2.47	2.96	3.45	3.94	4.44	4.93	5.42	5.92
NOM* 12	1.50	2.00	2.50	3.00	3.50	4.00	4.50	5.00	5.50	6.00
14	1.52	2.02	2.53	3.04	3.55	4.06	4.56	5.07	5.58	6.08
16	1.53	2.05	2.55	3.07	3.59	4.12	4.61	5.13	5.66	6.16
RBP** 18	1.55	2.07	2.58	3.11	3.64	4.18	4.67	5.20	5.74	6.24
RBP** 20	1.57	2.10	2.60	3.14	3.68	4.24	4.73	5.26	5.83	6.32
RBP** 22	1.59	2.12	2.63	3.18	3.73	4.30	4.78	5.33	5.91	6.40
24	1.60	2.14	2.65	3.22	3.78	4.36	4.84	5.40	-	-
26	1.62	2.17	2.68	3.25	3.82	4.42	-	-	-	-

*Nominal Pressure. The nominal in vitro device specifications do not take into account any lesion resistance.

**Rated Burst Pressure - Do Not Exceed (1.5 - 4.0 mm: 22 atm, 4.5 - 5.0 mm: 20 atm, 5.5 - 6.0 mm: 18 atm)

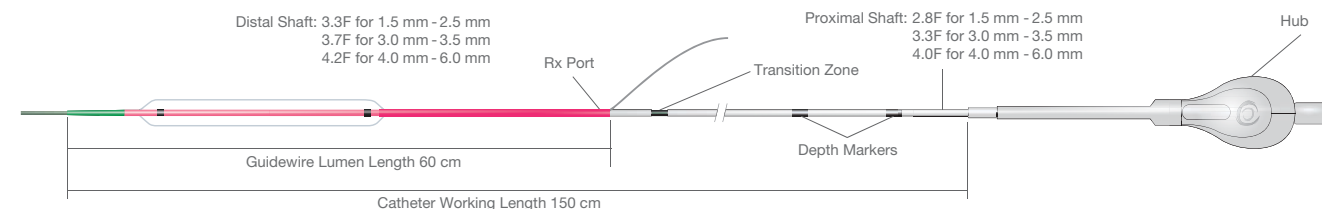

For more information please visit our website at www.OrbusNeich.com or contact us:

Corporate Headquarters

Units 303 & 305, 3/F, Building 20E
Hong Kong Science Park
Shatin, N.T., Hong Kong, China
Phone +852.2802.2288
Fax +852.2507.3532

EMEA Regional Headquarters

Drs. W. van Royenstraat 5
3871 AN Hoevelaken
The Netherlands

Phone +31.33.254.1150

Toll Free Phone 00800.0254.1150*

* Only for Belgium, Denmark, France, Germany, Ireland, Netherlands, Norway, Sweden and UK

©2025 OrbusNeich Medical Group Holdings Limited or its affiliates. All rights reserved.
The JADE product is now approved in the USA and distributed with USA specific catalog numbers.


G-70-0546 Rev07

JADE and OrbusNeich are registered trademarks of OrbusNeich Medical Group Holdings Limited or its affiliates.
All data and photos on file. Illustrations are not to be considered as engineering drawings or photos.

180 mm & 240 mm
balloon length
NOW
AVAILABLE!

Optimal Combination of
Robustness and Superb Deliverability

Strength and precision without compromising deliverability

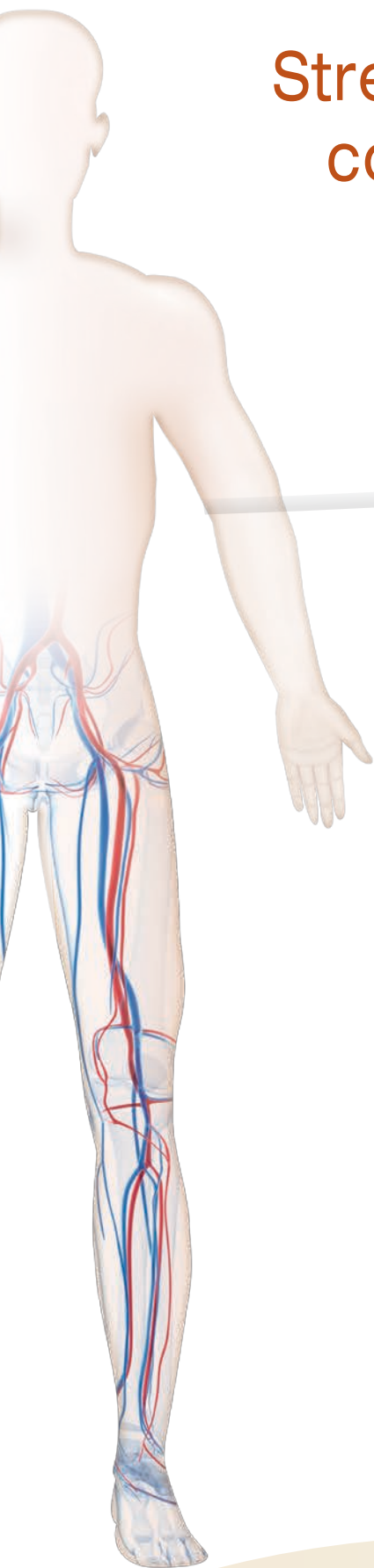

Z-tip technology for managing challenging lesions

High burst pressure **22 RBP** for opening resistant lesions

Extended guidewire tracking lumen facilitates extra pushability to reach very distal lesions.


Technical Specifications

Catheter Type	Rapid Exchange
Minimum Sheath Size	4F (Ø 1.5 mm - 3.5 mm); 5F (Ø 4.0 mm - 5.5 mm, Ø 6.0 mm for lengths 15 to 80 mm); 6F (Ø 6.0 mm for lengths 120 to 240 mm)
Guidewire Compatibility	0.014"
Catheter Working Length	150 cm
Proximal Shaft Diameter	2.8F (Ø 1.5 mm - 2.5 mm); 3.3F (Ø 3.0 mm - 3.5 mm); 4.0F (Ø 4.0 mm - 6.0 mm)
Distal Shaft Diameter	3.3F (Ø 1.5 mm - 2.5 mm); 3.7F (Ø 3.0 mm - 3.5 mm); 4.2F (Ø 4.0 mm - 6.0 mm)
Balloon Material	Nylon / Pebax
Compliance	Controlled compliance
Leading Tip Length	3.0 mm (Ø 1.5 mm - 3.0 mm); 4.0 mm (Ø 3.5 mm - 4.5 mm); 5.0 mm (Ø 5.0 mm - 6.0 mm)
Marker Bands	2 marker bands (Ø 15 mm - 120 mm); 4 marker bands (Ø 180 mm & 240 mm)
Coating	Hydrophilic (tip & distal shaft); Hydrophobic (guidewire lumen & balloon)
Crossing Profile	0.039" (Ø 3.0 mm)
Nominal Pressure	12 atm
Rated Burst Pressure	22 atm (Ø 1.5 mm - 4.0 mm); 20 atm (Ø 4.5 mm - 5.0 mm); 18 atm (Ø 5.5 mm - 6.0 mm)

True Non-compliance

Competitor Compliances Comparison for BTK* (2.0 mm)

Balloon Diameter (mm)

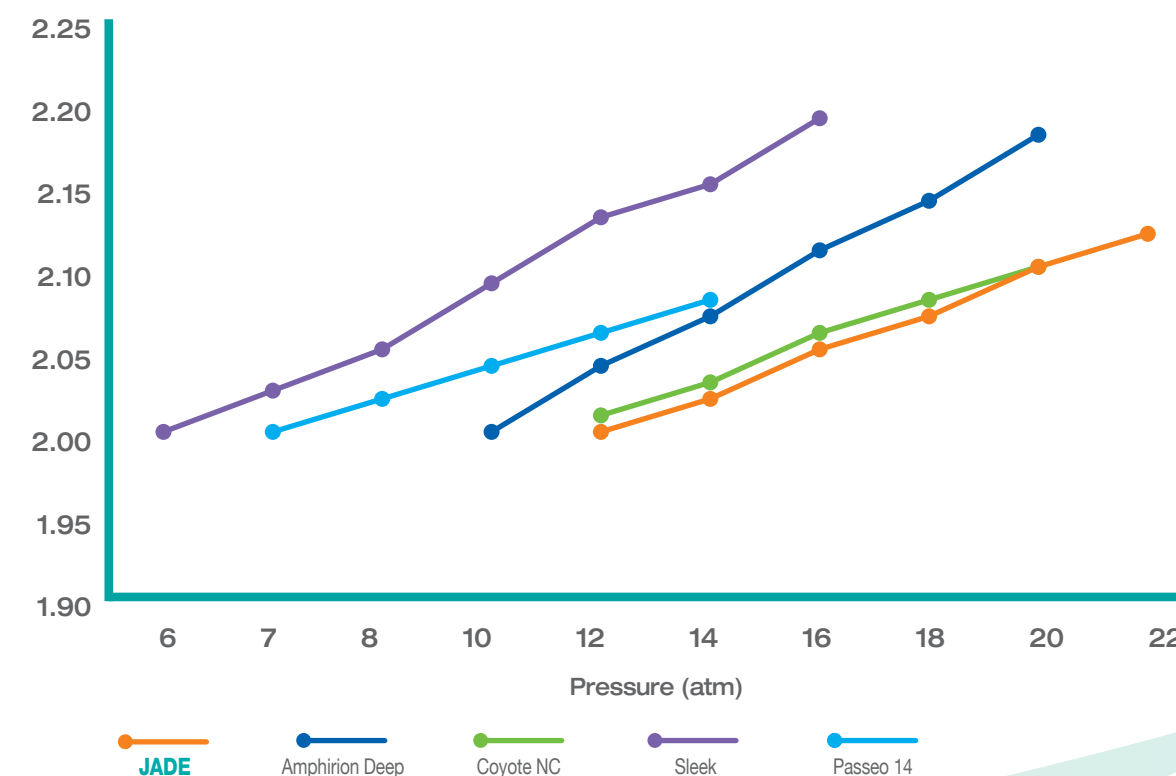

* Based on respective competitor balloons growth rate from NOM to RBP